

FLX Floorline Bed

Refined Floorline Care



FLX Floorline Bed



NEW
Customisable
Entry & Exit
Height Feature
 for improved safety
 during transfers

Shown with optional Como Curved head and foot boards



Included Features

- 1 Respite Ready 275kg Safe Working Load
- 2 UltraLow 100mm to 750mm Height Range
- 3 Single and King Single Sizes
- 4 Four Section Wipe Clean Mattress Platform
- 5 Adjustable Powered Backrest
- 6 Adjustable Powered Legrest
- 7 Anti-Shear Auto-Regression Backrest
- 8 Trendelenburg / Reverse Trendelenburg
- 9 Backlit Handset
- 10 Lockout Controls on Nurse Control Panel
- 11 Gentle Halo Underbed Lighting
- 12 Wall Buffer Bar
- 13 Integrated Bed Extension from 100-200mm
- 14 Side and Foot-End Mattress Retainers
- 15 One-Touch Foot-End Brakes with On / Off Highlights
- 16 Stylish, Custom and Bespoke Bed Ends and Side Panels
- 17 One Button Chair Positioning
- 18 One Button Entry / Exit Height
- 19 Adjustable Height Settings
- 20 Disassembles into two sections, offering improved storage and portability

Best Practice Care

Safety First

- **Custom Entry/Exit Height Function:** Set the exact safe transfer height for each resident. Effective across a large range of user heights.
- **Ultra Low Height:** Market leading 100mm low height whilst maintaining full ground clearance.
- **Attendant Control Panel:** Allows caregivers functional control and the ability to lock out certain bed functions from being operated via the handset.
- **Underbed Clearance:** Full clearance under the bed minimises the risk of cord impingement and protects flooring from damage.
- **SmartDescent:** Bed lowering speed is reduced once the bed reaches 300mm maximizing carer potential to remove obstacles as the bed lowers to the floor.
- **Underbed Lighting:** Creates a gentle halo around the bed to reduce falls risk and improve carer safety.

Clinical Considerations

- **Backrest Auto-Regression:** Reduces shear forces on the spine and sacral regions. Minimizes abdominal compression when back is raised supporting cardio pulmonary and respiratory function.
- **One Touch Chair Positioning:** Improves social interaction, assists respiratory and cardio pulmonary function and creates a comfortable position for sitting, dining and transfers.
- **Integrated Bed Extension:** 100-200mm length extension is easily operated and ensures a large range of user heights will fit comfortably without compromise.
- **Trendelenburg (and Reverse):** Augments patient positioning, comfort and assists with bed mobility*.

*Carers should seek advice from a trained Allied Health professional for correct manual handling techniques.

Infinite Positioning



HI LO



BACK RAISE



LEG RAISE



SEATED



TRENDELENBURG



Low position, back and knee raised
Sublime Teak side panels (accessory)



Floorline position
Self help pole (accessory)



Floorline position
Split folding grab rails (accessory)

Interior Styling

Contemporary Designs

Choose from our contemporary standard bed ends and sides; made to order styles; or completely bespoke designs of your choosing to meet any interior design requirements.

The FLX bed offers a suite of aesthetic options from our subtle Como Radius bed ends (standard) and side panels, to custom shapes and finishes.

Standard Bed Ends – Como Radius



Sublime Teak



Milano Walnut

Standard Side Panels



Sublime Teak



Milano Walnut

Custom Designs



Sublime Teak

Tallo
Generous outer frame and recessed inner panel with void for a luxurious aesthetic and a softer silhouette. Timber or laminate.



Milano Walnut

Como Curved
Sweeping curves provide a gentle aesthetic to compliment a robust outer frame design. Timber or laminate.



Milano Walnut

Eyre
Elegant frame and recessed inner panel with curved void detail creates an aperture to the surrounding room. Timber or laminate.



Sublime Teak

Sophia
Luxurious upholstery and subtle lines provide a gentle aesthetic to compliment interior design accents. Timber or laminate.

Bespoke Finishes

Don't see a design that suits your needs? We can customise as required to incorporate your choice of material, colour, finish, size or shape and match any interior design specifications.

Performance and Durability

- **Quick Swap Actuators:** Allow for quick and easy inspection and maintenance, minimising downtime.
- **TiMotion Electronics:** Medical grade electronics from a leading global manufacturer.
- **Control Box:** Easily accessible beneath the back rest panel for efficient servicing.
- **Quick Connect Bed Ends:** Simple and quick to set-up upon bed installation. Easily removed for access or maintenance.
- **Wall Buffer:** Length adjustable to suit room configuration.
- **Width Extension Kit:** Ability to increase a single to a king single, ideal for resident changes in a facility.

- **EasyMove Transport System:** Facilitates safe no-lift transportation and movement of beds, as well as safe and easy access underneath the bed platform for routine inspection and maintenance.
- **Cable Conduit:** Allows other cables like air mattress power cords to be safely tucked away, preventing trip hazards and cord damage.
- **High Performance Twin Castors x 8:** Head end total locking (x2), foot end central locking (x4).
- **Respite Ready:** Generous 275kg SWL facilitates care of bariatric residents and augments standard bed assets to facilitate respite care.
- **Pull Apart Capability:** The bed can be disassembled into two sections, offering improved storage and portability.


Technical Specifications

Patient Surface – Single	2000 x 900mm
Patient Surface – King Single	2000 x 1060mm
Overall Length	2130 – 2330mm
Length Extension	100 / 200mm
Overall Width – Single	964mm
Overall Width – King Single	1124mm
Patient Surface Height Range	100 – 750mm
Safe Working Load	275kg
Backrest Range	0 – 62°
Knee Bend Range	0 – 40°
Underbed Clearance	Variable. Min 20mm
Mattress Platform	4 Sections
Product Weight – Single	105kg
Product Weight – King Single	115kg
Castors	8 x Twin 75mm

Electrical Parameters

Power Input	240V AC 50Hz, Max current 4000mA
Rate Load	Interruption 10%, Max 2 minutes / 18 minutes
Leakage Protection	Class II Type B
Waterproof Grade	IP54

Comprehensive user manual, troubleshooting guide, safety checklist and assembly instructions available.

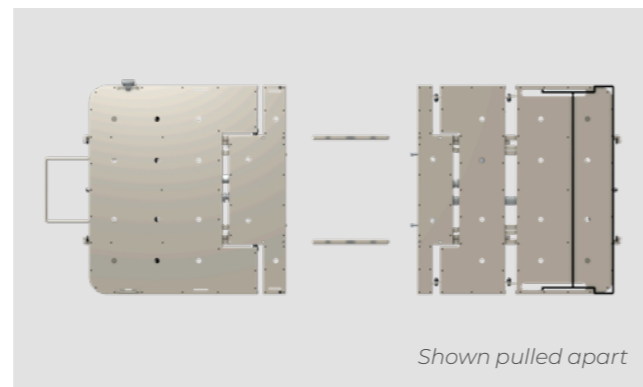

Beds

- BEB048600 FLX Floorline Bed – S
- BEB048610 FLX Floorline Bed – KS

Standard Side Panels (Sold Separately)

- BEB046265 Sublime Teak
- BEB046275 Milano Walnut

Standard Bed Ends (Sold Separately)

- BEB048700 Como Radius – Sublime Teak – S
- BEB048701 Como Radius – Milano Walnut – S
- BEB048710 Como Radius – Sublime Teak – KS
- BEB048711 Como Radius – Milano Walnut – KS


✓ FLX Bed is registered in the ARTC, AUST R ARTC

3 YEARS WARRANTY Electronics

10 YEARS WARRANTY Frame

1 YEAR WARRANTY Handset + Cables

Accessories

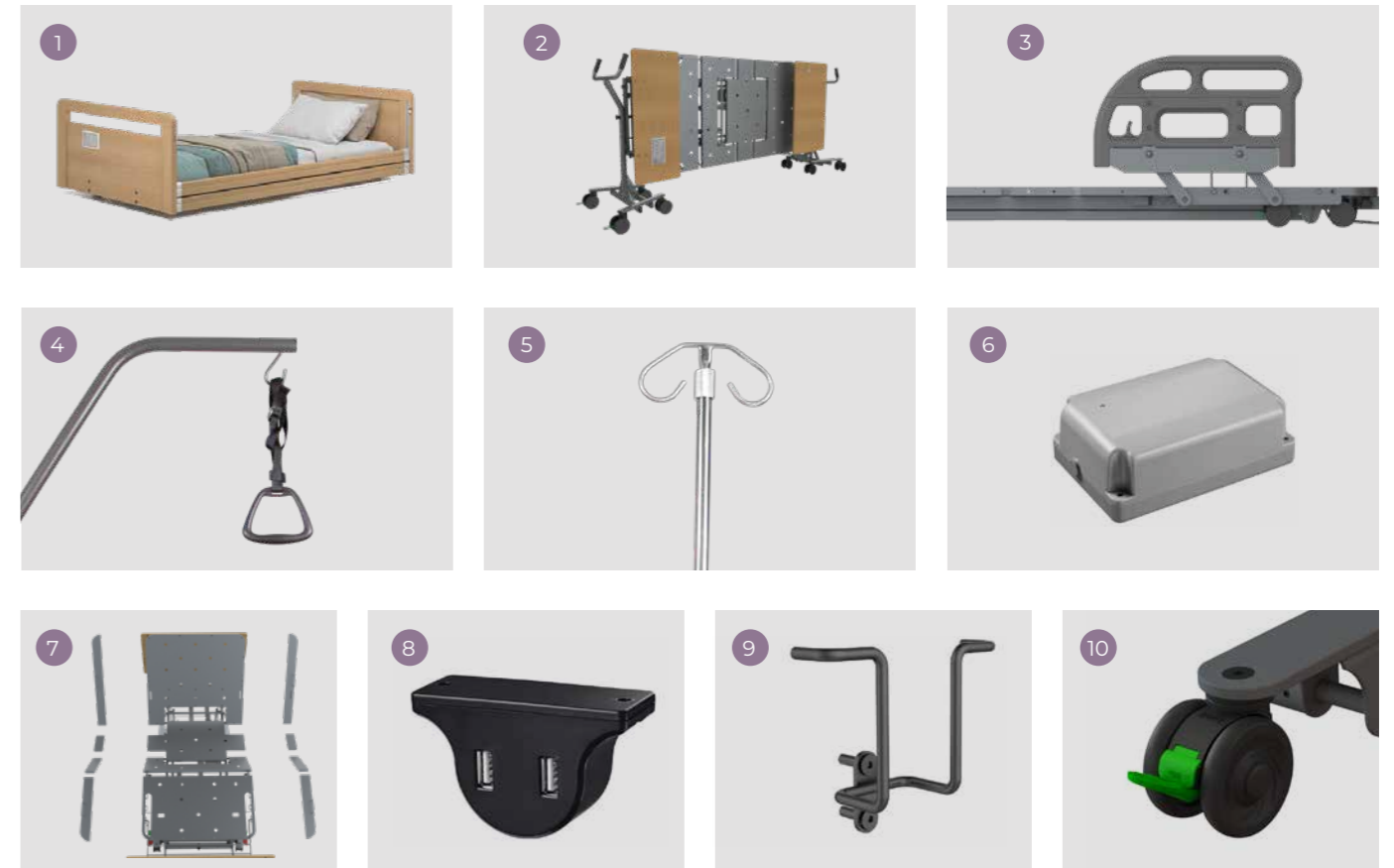

- 1. Full Length Rails**
BEB048785 - Sublime Teak
BEB048790 - Milano Walnut
End to end prevention of falls. Includes both left/right sides. Easily lifted by hand, collapse to lowest position beside the bed platform.
*Bed ends sold separately
- 2. EasyMove Transport System** BEA030055
Is an ideal transport solution designed to minimise manual handling. It allows effortless lifting and tilting of the bed to its transport position, while the large 100mm castors offer smooth and easy manoeuvring for one person. The EasyMove is also a perfect maintenance platform facilitating easy access to mechanical parts under the bed platform.
- 3. Split Folding Grab Rails (Pair)** BEB046525
Left/Right pair of moulded rails provides manoeuvring and repositioning assistance to a resident. Can be positioned at the head end or foot end. Folds down to lowest position beside the bed platform.
- 4. Self Help Pole** BEA829140
Self Help Pole Mounting Bracket BEA829145
Provides a firm support handle for improved bed mobility.
- 5. IV Pole** BEA029710
IV Pole Mounting Bracket BEA829145
- 6. Battery Backup** BEB046592
Prolong bed operation by up to 15 lifts during a power outage.
- 7. Width Extension Kit** BEB046500
Allows conversion of a Single FLX into a King Single.
- 8. USB Port Upgrade Kit** BEB046616
Add USB charging functionality to the side of the bed, perfect for charging mobile phones and tablets.
- 9. Handset Holder** BEA030030
Mounts discreetly behind the head board for easy access and improved room aesthetics.
- 10. Directional Castor** BEB046630
Replaces total locking castor at the head end for improved bed tracking during bed movement.



ABBR13 211124

 **aidacare**
FLX FLOORLINE BED

09 520 2962 | aidacare.co.nz