



SHELL CHAIR

For every occasion, easy to use in community or facility environments



aspire for...

Comfort

Support

Protection

Maximises Independence & Function

Lightweight Ergonomic Design

Customisable Postural Support

Genuine Dartex PU Fabric

Height Adjustable Seat Mount

aidacare.co.nz
09 520 2962

SHELL CHAIR

The Aspire Shell Chair is a comfortable mobile seating solution designed to maximise user access and social interaction within an Aged Care Facility or Healthcare environment.

The compact design allows access under tables for sharing meals, games and other activities. The lightweight frame with anterior tilt function for sit to stand transfers reduces the manual handling burden for caregivers.

TILT FUNCTIONALITY

Forward tilt assists with chair egress and standing transfers. Backwards tilt assists in pressure redistribution and patient stability during transfers.



COMFORT & SUPPORT

Ergonomic wing headrest and lumbar cushioning for improved seating posture

SAFETY LAP BELT

Supportive lap belt designed with a single pull strap for convenient adjustment. Can also be tucked out of the way when not in use

ADJUSTABLE SEAT HEIGHT

Seat height can be adjusted by up to 80mm allowing the chair to be used at a desk or dining table

ERGONOMIC LEG REST

The supportive leg rest automatically raises or lowers as the chair is tilted



ATTENDANT CARER PUSH HANDLE

Ergonomically designed, tilt trigger lever makes postural adjustment a simple task

REMOVABLE SEAT CUSHION

Facilitates acute pressure care options such as utilising a ROHO cushions

STORAGE POCKET


Conveniently carries everyday items such as clip boards, magazines or diaries

LIGHTWEIGHT AND MANOEUVRABLE

12.5" rear wheels with attendant foot brake, 5" front swivel castors together with a lightweight design, makes for easy and safe handling for caregivers and a smooth and comfortable ride for users

OPTIONAL ACCESSORIES

1. Pillow Headrest
CHP198930
2. Footplate
CHP198925
3. Retractable Lap Tray
CHP198920

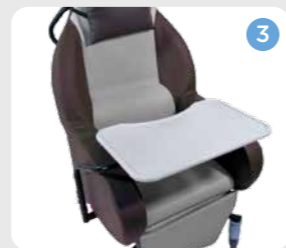

- DVA Contracted
- AS / NZS Approved
- Shell Chair is registered in the ARTG, AUST R ARTG 277232

Clinical Consideration

OPTIMISES USER INDEPENDENCE

- Anterior tilt functionality assists with chair egress and standing transfers, maximises a user's independence in a safe manner

ANTHROPOMETRIC CONSIDERATIONS

- Two sizes available, with height adjustment, to ensure correct fit and postural support

FATIGUE MANAGEMENT

- The 27° tilt allows users to relax whilst sitting unassisted

PRESSURE INJURY MANAGEMENT

- Dartex fabric is high stretch reducing the risk of shear and friction
- The tilting feature allows pressure redistribution

OEDEMA MANAGEMENT

- Raised leg positioning to assist with lower limb oedema
- Compatible with other, pressure relieving cushions

REDUCED MANUAL HANDLING BURDEN

- A lightweight design and small footprint to make manoeuvring the chair easier for caregivers
- The anterior tilt function means that a user can be more independent in sit-to-stand transfers, reducing the need for more than one carer to assist
- As a result, the need for manual handling equipment is reduced



1 YEAR WARRANTY
Fabric


2 YEAR WARRANTY
Frame

CODE	DESCRIPTION
CHP198910	Shell Chair - Medium
CHP198915	Shell Chair - Large

GENUINE DARTEX® PU FABRIC

Well respected in the healthcare industry, Dartex® PU fabric is the perfect patient surface solution for treatment and management of pressure injuries through it's 5 key characteristics (Figure 1).

FIGURE 1

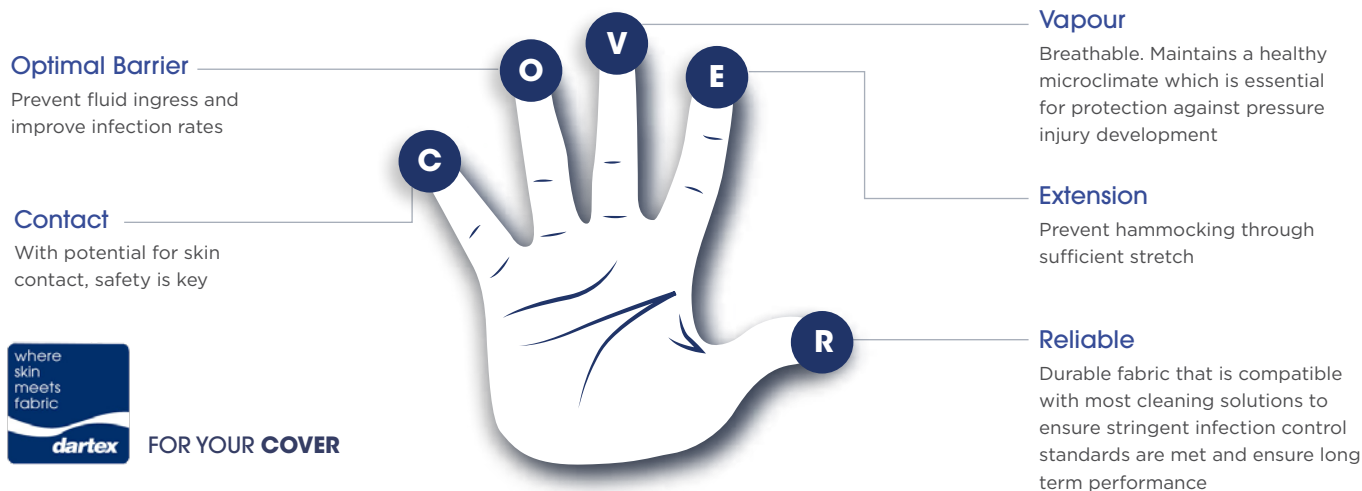

PATIENT & CARER BENEFITS



Designed specifically for pressure care needs



Removable anti-bacterial/anti-microbial cover for infection control management



Waterproof yet breathable for improved skin microclimate

TECHNICAL SPECIFICATIONS

FEATURES	M	L
Seat Height (mm):	450 - 530	450 - 530
Seat Depth (mm):	440	470
Seat Width (mm):	450	510
Armrest Height (mm):	230	230
Overall Width (mm):	700	720
Overall Length (mm):	890	890
Overall Weight:	30kg	30kg
Tilt Angle:	-3° to +27°	-3° to +27°
SWL:	135kg	135kg

* Chair measurements may vary +/- 10mm due to upholstery tolerances
* Seat width is taken between the thickest sections of armrest padding

EXTERIOR SURFACE

Brown (leather-style) healthcare fabric has waterproof, anti-fungal, anti-bacterial and fire retardant properties.

PATIENT SURFACE

Dartex Cool Grey

- High stretch fabric for pressure reduction
- Vapour permeable (breathable)
- Anti-bacterial
- Anti-fungal
- Fire retardant

CUSHION SURFACE

- Removable seat cushion facilitates pressure care options
- Compatible with after-market seating such as ROHO cushions