

# FL250 Floorline Bed

Low Height, High Style





# FL250 Floorline Bed



SWL  
250KG

\*Shown with Tallo Bed Ends (Custom Made)

## FL250 Features




1. 250 SWL
2. 110mm Low Height / 740mm Max Height
3. Single and King Single Sizes
4. Four Section Wipe Clean Mattress Platform
5. Adjustable Powered Back Rest
6. Adjustable Powered Leg Rest
7. Anti-Shear Auto-Regression Backrest and Legrest
8. Trendelenburg / Reverse Trendelenburg
9. Handset - Backlit Attendant Control Panel - Lockable
10. Gentle Halo Underbed Lighting
11. Intelligent Wall Protection Buffer
12. Integrated Bed Extension from 100mm and 200mm
13. Side and Foot-End Mattress Retainers
14. In-Line Central Locking Brakes
15. Stylish, Custom and Bespoke Bed Ends and Side Panels
16. One Button Chair Positioning
17. One Button Entry / Exit Heights
18. Adjustable Height Settings



# Best Practice Care

## Safety First

- **Low Height:** With its impressive low height of 110mm, injury to the resident can be significantly reduced or avoided in the event of a fall from bed, therefore reducing the need for side-rails.
- **Attendant Control Panel:** Allows caregivers functional control and the ability to lock out certain bed functions from being operated via the handset.
- **Underbed Clearance:** Market leading 110mm low height whilst maintaining full ground clearance.
- **ExitSafe:** By pushing the ExitSafe button the FL250 moves to a flat and ergonomic position at the optimum height for resident ingress / exit.
- **BrakeFirst:** Brakes automatically lock when bed reaches its lowest height. Brakes are located within the footprint of the bed, minimising any trip hazard for the caregiver and resident alike. Brakes are easily accessible at all heights.
- **SmartDescent:** Bed lowering speed is reduced once the bed reaches 300mm maximising carer potential to remove obstacles as the bed lowers to the floor.
- **Underbed Lighting:** Creates a gentle halo around the bed to reduce falls risk and improve carer safety.

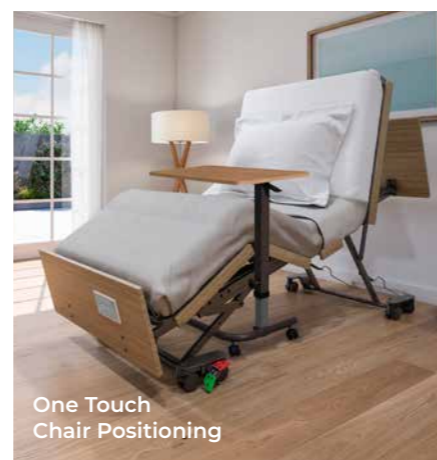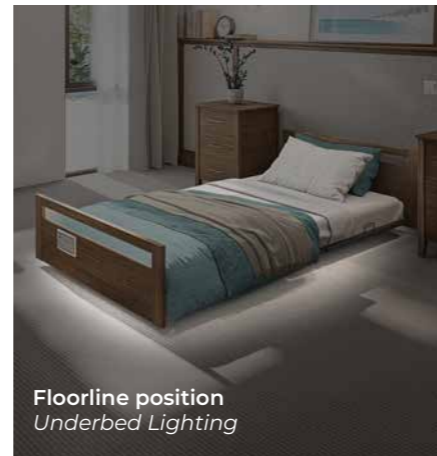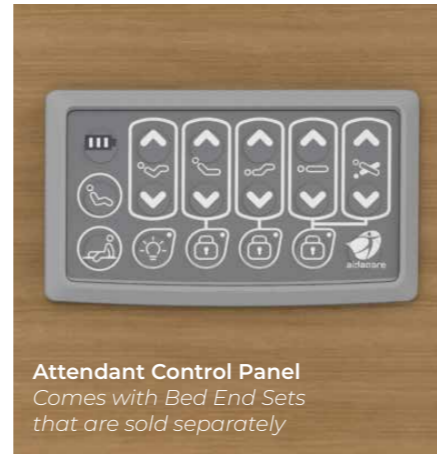

## Clinical Considerations

- **Sleeping Low is Clinically Proven:** International clinical studies have shown the benefits of sleeping low and eliminating side rails to resident safety, welfare and comfort, in particular by lowering fall rates and fall injury severity.
- **Auto-Regression:** AutoRegression helps to reduce sacral pressure, shear forces and abdominal compression during bed profiling.
- **CPR Position:** The low height and flat surface of the ExitSafe position also facilitates the administration of CPR in the event of an emergency.
- **One Touch Chair Positioning:** Improves social interaction, assists respiratory and cardio pulmonary function and creates a comfortable position for sitting, dining and transfers.
- **Integrated Bed Extension:** 100-200mm length extension is easily operated and ensures a large range of user heights will fit comfortably without compromise.
- **Trendelenburg (and Reverse):** Augments patient positioning, comfort and assists with bed mobility\*.

\*Carers should seek advice from a trained Allied Health professional for correct manual handling techniques.

## Infinite Positioning



# Interior Styling

## CONTEMPORARY DESIGNS

Choose from our contemporary standard bed ends sides and rails; made to order styles; or completely bespoke designs of your choosing to meet any interior design requirements.

The FL250 bed offers a suite of aesthetic options from our subtle Como Radius bed ends (standard) and side panels, to custom shapes and finishes.

## Standard Bed Ends – Como Radius



## Standard Side Panels



## Custom Bed End Designs

Range of options to achieve desired aesthetic, silhouettes, elegance and luxurious styles. Available in timber or laminate finishes as well as various colours



**Bespoke Finishes:** Customise material, colour, finish, shape or size to match any interior design specifications.

## Side Rails

The FL250 offers a range of contemporary side rail options to meet patient safety requirements needs without compromising on style and comfort.

**TELESCOPIC RAILS:** One-handed operation to lift and lower rails. Rails discretely lower behind laminate boards when not in use. 1/3, 2/3 or full length configurations can be mounted individually or together. Available in Sublime Teak or Milano Walnut.

*The minimum bed height changes to 275mm when Telescopic Rails installed*

**FULL LENGTH RAILS:** Rails neatly lower when not in use to create a single uninterrupted timber-look surround to the bed. Available in Sublime Teak or Milano Walnut.

*The minimum bed height remains unchanged when Full Length Rails installed*



Telescopic rails

# Performance and Durability

- **Quick Swap Actuators**  
Allow for quick and easy inspection and maintenance, minimising downtime.
- **TiMotion Electronics**  
Medical grade electronics from a leading global manufacturer.
- **Control Box**  
Easily accessible beneath the back rest panel for efficient servicing.
- **Quick Connect Bed Ends**  
Simple and quick to set-up upon bed installation. Easily removed for access or maintenance.
- **Wall Buffer**  
Length adjustable to suit room configuration.
- **Cable Conduit**  
Allows other cables like air mattress power cords to be safely tucked away, preventing trip hazards and cord damage.
- **EasyMove Transport System**  
Facilitates safe no-lift transportation and movement of beds, as well as safe and easy access underneath the bed platform for routine inspection and maintenance.
- **Width Extension Kit**  
Ability to increase a single to a king single, ideal for resident changes in a facility.
- **High Performance Twin Castors x 8**  
Head end directional locking (x2), foot end central locking (x4)
- **Respite Ready**  
Generous 250kg SWL facilitates care of bariatric residents and augments standard bed assets to facilitate respite care.

# Specifications

## Technical Specifications

Patient Surface – Single	2000 x 900mm
Patient Surface – King Single	2000 x 1060mm
Overall Length	2088 – 2288mm
Length Extension	100 / 200mm
Overall Width – Single	900mm
Overall Width – King Single	1060mm
Patient Surface Height Range	110 – 740mm
Safe Working Load	250kg
Backrest Range	0 – 60°
Knee Bend Range	0 – 35°
Underbed Clearance	Variable. Min 30mm
Mattress Platform	4 Sections
Product Weight – Single	105kg
Product Weight – King Single	115kg
Castors	8 x Twin 75mm

## Electrical Parameters

Power Input	230V AC 50Hz, Max current 1A
Rate Load	Interruption 10%, Max 2 min/18 min
Leakage Protection	Class I Type B
Waterproof Grade	IP54

## FL250 Floorline Beds

BEB046100 Single	BEB046150 King Single
------------------	-----------------------

## Standard Side Panel Sets (Sold Separately)

BEB046240 Sublime Teak	BEB046250 Milano Walnut
------------------------	-------------------------

## Standard Como Radius Bed End Sets (Sold Separately)

BEB048700 Sublime Teak Single
BEB048701 Milano Walnut Single
BEB048710 Sublime Teak King Single
BEB048711 Milano Walnut King Single

Comprehensive user manual, troubleshooting guide, safety checklist, assembly instructions available.

**10 YEARS WARRANTY**  
Frame

**3 YEARS WARRANTY**  
Electronics


**1 YEAR WARRANTY**  
Handset + Cables

FL250 Bed is registered in the ARTG, AUSTRALIA

# Accessories



- 1. Bed Ends for Telescopic Rails - Como Radius**  
BEB048830 - Sublime Teak - Single  
BEB048835 - Milano Walnut - Single  
BEB048840 - Sublime Teak - King Single  
BEB048845 - Milano Walnut - King Single
- 2. Bed Ends for Full Length Rails - Como Radius**  
BEB048750 - Sublime Teak - Single  
BEB048755 - Milano Walnut - Single  
BEB048770 - Sublime Teak - King Single  
BEB048775 - Milano Walnut - King Single
- 3. Full Length Rails**  
BEB048785 - Sublime Teak  
BEB048790 - Milano Walnut  
End to end prevention of falls. Includes both left/right sides. Easily lifted by hand, collapse to lowest position beside the bed platform.  
\*Bed ends sold separately
- 4. Telescopic Rails (Full Length: 1/3 + 2/3)**  
BEB048800 - Sublime Teak  
BEB048805 - Milano Walnut  
Adjustable height rails to assist in the prevention of falls as well as manoeuvring. Includes both left/right side rails, and bed height limiting handset. Various configurations include 1/3, 2/3 individually or together for full length. \*Bed ends sold separately
- 5. Telescopic Rails (Part Length + Decorative Trims)**  
BEB048810 - 1/3 rails + 2/3 DT - Sublime Teak  
BEB048815 - 1/3 rails + 2/3 DT - Milano Walnut  
BEB048820 - 2/3 rails + 1/3 DT - Sublime Teak  
BEB048825 - 2/3 rails + 1/3 DT - Milano Walnut
- 6. Bed Rail Protector - Full Length Side Rails**  
BEA010778 - Mesh

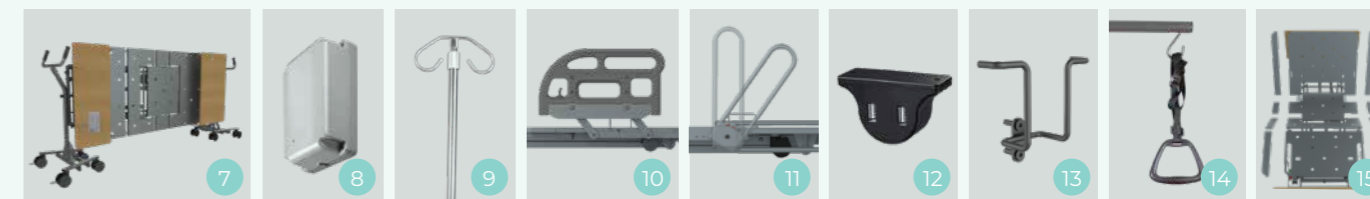

- 7. EasyMove Transport System BEA030055**  
Is an ideal transport solution designed to minimise manual handling. It allows effortless lifting and tilting of the bed to its transport position, while the large 100mm castors offer smooth and easy manoeuvring for one person. The EasyMove is also a perfect maintenance platform facilitating easy access to mechanical parts under the bed platform.
- 8. Battery Backup BEB046592**  
Prolong bed operation by up to 15 lifts during a power outage.
- 9. IV Pole BEB046530**
- 10. Split Folding Grab Rails (Pair) BEB046525**  
Left/Right pair of moulded rails provides manoeuvring and repositioning assistance to a resident. Can be positioned at the head end or foot end. Folds down to lowest position beside the bed platform.
- 11. Folding Assist Bar (Single) BEB046520**  
Provides manoeuvring assistance to a resident. Folds down when not in use, out of the way of transfers.
- 12. USB Port Upgrade Kit BEB046616**  
Add USB charging functionality to the side of the bed, perfect for charging mobile phones and tablets.
- 13. Handset Holder BEA030030**  
Mounts discreetly behind the head board for easy access and improved room aesthetics.
- 14. Self Help Pole BEB046540**  
Provides a firm support handle for improved bed mobility.
- 15. King Single Conversion Kit BEB046500**  
Allows conversion of a Single FL250 into a King Single.





ABBR07 211124

 **aidacare**  
FL250 FLOORLINE BED

09 520 2962 | [aidacare.co.nz](http://aidacare.co.nz)